




Proven Products for Better Dentistry




PRODUCT BROCHURE

# Select from a Rational Shade Line Up


## Enamel

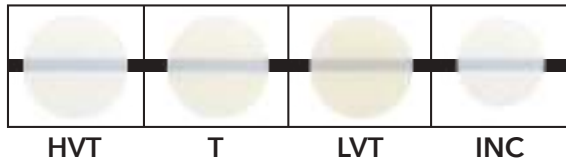


## Bleach White

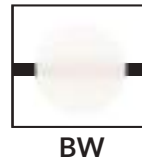


## Universal

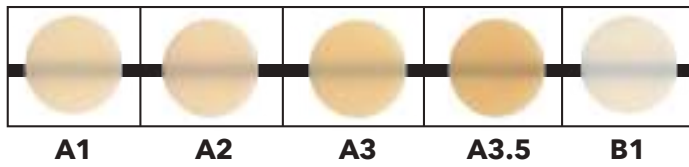


## Dentin

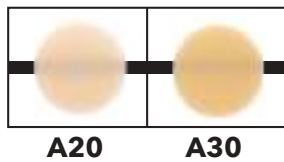



## Opaquer

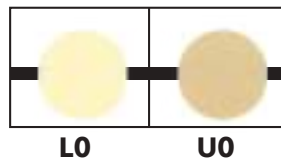

## Bulkfill

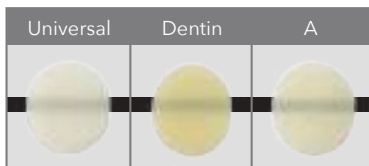

## Gingiva

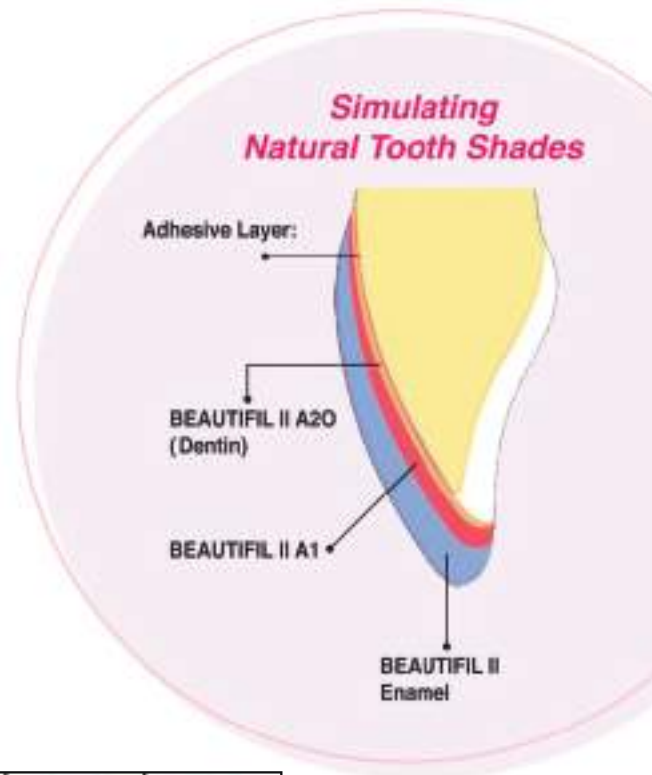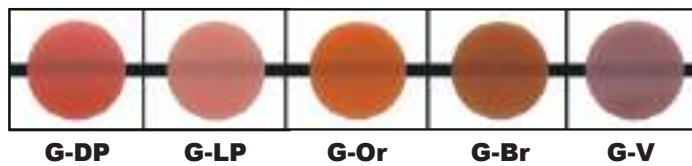

### GUIDE FOR CLASS III RESTORATION WITH ENAMEL SHADES

Tooth Preparation	Starburst Bevel Technique	Lingual Shelf, Enamel	Caulking the Window Technique	Dentin Replacement	Enamel Layer	Final Restoration

### GUIDE FOR CREATING ARTISTIC COMPOSITE VENEER WITH ENAMEL SHADES

Tooth Preparation	Lingual Enamel	Dentin Layer	Dentin Layer, Mamelons, Incisal Halo	Cervical Enamel	Incisal Enamel	Final Restoration

## Beautiful II LS



- Bioactive, Nano-Hybrid Composite.
- Great strength and high wear resistance.
- Polishes in an instant.
- Tooth-like aesthetics with natural fluorescence and chameleon effect.
- Low volumetric shrinkage (0.85%) and shrinkage stress (2.72 MPa).
- Available as 4 gm syringes, and 0.25 capsules.



## Beautiful II



- Bioactive, Nano-Hybrid Composite.
- Light transmission and diffusion properties of natural teeth, and great chameleon effect.
- Reliable strength and longevity.
- Exceptional radiopacity.
- Available as 4.5 gm syringes, and 0.25 gm capsules.



## Beautiful Flow Plus & Plus X



- Bioactive flowable/injectable nano-hybrid composite. Indicated for all classes-including occlusal and cusp.
- Improved handling and effortless polishing and maintains a long-lasting shine.
- Superior adaptation ensures a tight marginal seal, decreasing voids.
- F00 – Zero Flow: Excellent shape retention, ideal for build-up: The flowability of F00 is so low that the paste is not deformed during layering.
- F03 Low Flow: Moderate shape retention and viscosity plus perfect adaptation: The self-levelling F03 paste thoroughly wets the cavity walls and is therefore ideal for filling.
- Available as 2.2 gm syringes.

Beautiful Flow Plus



Beautiful Flow Plus X



## Beautiful-Bulk Restorative



- Packable/paste Bulk fill composite designed to deliver aesthetic results, low shrinkage stress and complete polymerization at 4 mm depth of cure.
- High fill ratio (87.0 wt% 74.5 vol%) reduces polymerization shrinkage and shrinkage stress while increasing compressive and flexural strength.
- Available as 4.5 gm syringes, and 0.25 gm capsules.



## Beautiful-Bulk Flowable



- Flowable Bulkfill composite designed to deliver aesthetic results, low shrinkage stress and complete polymerization at 4 mm depth of cure.
- Self-leveling allows easy adaptation.
- Optimal shade match (chameleon effect).
- Available as 2.5 gm syringes, and 0.23 gm capsules.



## Beautiful II Enamel



- For creating polychromatic restorations with morphology and aesthetics of natural teeth.
- Can be layered and blended to deliver results that will mimic natural dentition.
- Achieve highly aesthetic direct restorations with minimal shade selections.
- Available as 2.5 gm syringes.



## Beautiful II Gingiva



- Indicated in cases of gingival recession, exposed abutments and crown margins, root erosion and missing papillae
- Shades can be layered and combined for aesthetic solutions without the need for surgery
- Perfect for Class V restorations



## Beautiful Opaquer



- Mask discolored or stained teeth, metal color & pulp capping agent
- Light-curing, fluoride releasing restorative material
- Available in the two shades, Universal Opaque & Light Opaque




## FL-Bond II




- Self-etching fluoride releasing two-step adhesive system which features an excellent bond to both enamel and dentine and provides a secure marginal seal.
- High filler content.
- Adhesive bonding of any composite to tooth structure.
- Used also for Pre-treatment of the tooth structure prior to the application of resin cements.
- Available as a kit (5 ml primer + 5 ml adhesive).

Low viscosity (other bonding agent)      Moderate viscosity (FL-Bond II)

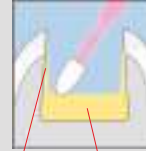

Extremely thin bonding layer      Collection of liquid at cavity floor

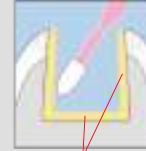

Equal bonding layer thickness




## BeautiBond Xtreme




- The powerful all-in-one universal adhesive for all etching techniques and a wide variety of restorative materials.
- Excellent bond strength, low technique sensitivity and convenient application in one step.
- Users no longer require an activator or enhancer for certain restorative materials.
- It does not matter whether tooth surfaces are wet or dry or what etching technique is used.
- The nozzle has been specially designed to deliver very small drops (approx. 260 per bottle), ensuring economical dispensing with only minimal loss of material.
- The adhesive film is only approx. 5  $\mu$ m thick, which makes this filler- and HEMA-free adhesive the system of choice, particularly for minimally invasive dentistry and highly aesthetic restorations.
- Available as 5 ml bottle.



Self-etch

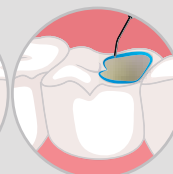

Selective-etch

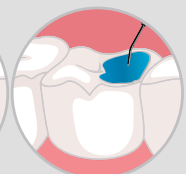

Total-etch



## Beauticem SA

- Durable bond to a wide range of substrates
- No need for separate phosphoric acid etching, priming or bonding of enamel and dentine
- Outstanding physical properties
- High translucency for aesthetically challenging cases
- Low film thickness (11.8 µm)
- Effective fluoride release for sustained anticariogenic and antibacterial properties



## Beauticem Veneer

- Intuitive shade selection
- Shade stability thanks to amine-free composition
- Outstanding ease of use
- Fluoride release and recharge
- High radiopacity
- Try-in Paste in the same shades
- Optimal shade matching of Try-in Paste and cured cement



## BeautiLink SA

- Self-Adhesive Resin Cement Specifically for Zirconia Restorations with Bioactive Giomer Technology
- High Bond Strength with No Primer Needed
- Automix Syringe & Ideal Handling
- Dual Cure-Fast & Easy Clean Up



## CX-Smart

- Chemically bonds to the tooth structure for excellent marginal seal and sharp set with remarkable moisture tolerance.
- No etch or bond for a quick, easy and simple application- ideal for paediatric, geriatric and special needs patients.
- Resin-free, suitable for patients with resin allergy.
- Optimal thixotropy and even film thickness with durable physical properties, high acid resistance and low oral solubility.
- HY-additive ensures virtually no post-operative sensitivity & favourable pulp response.
- Cost-effective and efficient.



## CX-Plus

An enhanced glass ionomer luting cement with prolonged working time.

- High resistance to compression 149 MPa compared with competitive cements
- Low film thickness, 15 microns optimizing the marginal integrity
- Radiopacity 1,55 mm
- Low level of acid erosion 0,02 mm
- High fluoride release
- Capsules with 0,53 g
- Working time 4 minutes and 10 seconds from start of mixing
- 3 years of shelf life



## Zirconomer Improved

- Remarkable edge strength and excellent marginal adaptation
- Self-adhesive with tooth-like co-efficient of thermal expansion
- Sustained fluoride protection, especially in cases with high caries risk
- Excellent mixing & handling qualities
- Higher translucency for a closer match to natural tooth
- Durability like silver amalgam without the hazard of mercury
- Excellent resistance to abrasion & erosion



## Glasionomer FX Ultra

- Remarkable Translucency For An Easier Shade Match And Better Aesthetics
- Superior mechanical properties for durable posterior restorations
- Good wettability of powder and liquid assures an easy and fast mix
- Superior Handling qualities with adequate working time and a sharp setting reaction
- Low solubility and acid erosion
- Significant fluoride release and recharge
- Easy to pack and adapt
- Available in shades A2, A3, A3.5, B2



## ZiLMaster

- Designed to contour, adjust, finish, and polish modern ceramic CAD/CAM materials, including zirconia, lithium silicate, and disilicate.
- Deliver long-lasting performance with a precise cut smooth finish, and tooth-like luster.
- High gloss in three polishing steps.
- No polishing paste is needed.



## Compomaster

- CompoMaster Coarse polishers have been designed for contouring and prepolishing composite surfaces in one step.
- CompoMaster instruments subsequently polish the restorations.



## Dura White

- Precise concentricity
- Low vibration grinding
- Consistent finishing performance
- Economical
- High edge stability
- Long-lasting service life
- Wide range of shape



## OneGloss

- One-step finisher and polisher that provides the ideal finish for all types of composite restorations.
- Simply by altering the contact pressure as you go from finishing to polishing - without instrument change and polishing paste.
- Available in 4 different shapes (cup - flame /midi point disc - mini point).
- Use intermittent water spray when applying high pressure or polishing enamel.



## Direct Dia Paste

- Diamond impregnated polishing paste for intra and extraoral high gloss polish.
- Paste contains 20% diamond particles in grain sizes of 2-4  $\mu\text{m}$ .
- It creates high gloss on natural enamel and all restorative materials (porcelain, direct and indirect composites, metals etc.).
- Available as 3 gm syringe.



## Super-Snap & Super-Snap X-Treme

- Super-Snap is designed for easy, fast and safe contouring, finishing and polishing of all microfilled and hybrid composites.
- Elastic silicone shank mount for easy placement of the disk onto the mandrel
- 3D semispherical shaped grits allow space to discharge ground debris while polishing, helping to reduce clogging and secondary scratches.
- Very flexible, ultra-thin disks allow easy access to inter-proximal areas.
- Available as assorted kits and single refills.



## Super Snap Buff Disk

- Polishing paste adheres to the felt coating
- Clean and economical working with polishing paste
- No metal centre



## Beautisealant



- Ideal handling and viscosity: eliminate voids with bubble-free consistency
- Radiopacity equal to dentine: 0.92 mm Al
- Unique cariostatic properties of S-PRG fillers
- High shear bond strength without phosphoric acid etch or rinse required
- Four easy steps; 30 second application
- Bisphenol A- and HEMA-free



## Pro-Care Gel



- Ion releasing treatment paste - PRG Pro-Care Gel α
- Treatment paste with S-PRG filler has appeared in preventive care



## PRG Barrier Coat



- Blocks and protects from external irritants for up to 6 months
- Fluoride release and recharge ability
- Excellent self-adhesive properties allows thin coating on flat and slanted surfaces (15um) - aesthetic results
- Acid neutralization and anti-plaque effect
- Simple and quick procedure
- HEMA-free, Acetone-free, and Alcohol-free
- Indicated for a variety of applications - orthodontic brackets, hard to brush areas, exposed root surface, newly erupted molars, white spots and areas of hypersensitivity



### SHOFU DENTAL MIDDLE EAST AND AFRICA

BUILDING 6EB, OFFICE 535, DUBAI AIRPORT FREEZONE, DUBAI, UAE, PO BOX: 371451  
TEL: +971 454 836 56 | EMAIL: info@shofudental.ae | www.shofu.com/global